

PERIODONTICS

Sharpen free LM Sharp Diamond™ coated instruments

Sickle scalers | Universal curettes | Dual Graceys |
Special instruments | Gracey curettes |
LM-ErgoMix™ implant instruments |
Periodontal hand instrument sets

DIAMOND FILES

Diamond coated files are designed for cleaning deep concave root surfaces in anatomically difficult areas.


Can be used with multidirectional strokes.



Concavity Diamond File SV 5-6

- Diamond coated round curved toe
- For acute furcation lesions


- LMErgoSense®** LM 255-256D ES
- LMErgoSense®** LM 255-256D EST
- LMErgoMax™** LM 255-256D XSI
- LMErgoNorm™** LM 255-256D SI

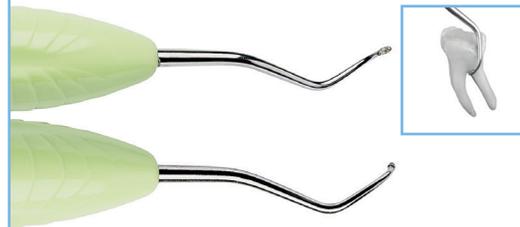

Furcation Diamond File SV 7-8

- Diamond coated egg-shaped toe
- For deep furcation lesions


- LMErgoSense®** LM 257-258D ES
- LMErgoSense®** LM 257-258D EST
- LMErgoMax™** LM 257-258D XSI
- LMErgoNorm™** LM 257-258D SI

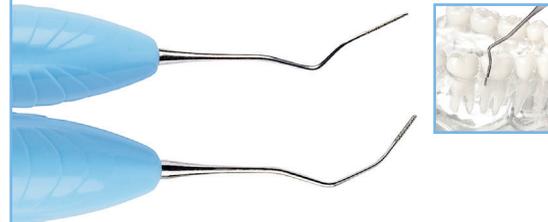

Diamond File mes-dist

- Extremely delicate diamond coated instrument for definitive scaling of the root surface
- Easy access with a fine, flexible blade in mesial-distal line angles and developmental grooves

- LMErgoSense®** LM 263-264D ES
- LMErgoSense®** LM 263-264D EST
- LMErgoMax™** LM 263-264D XSI

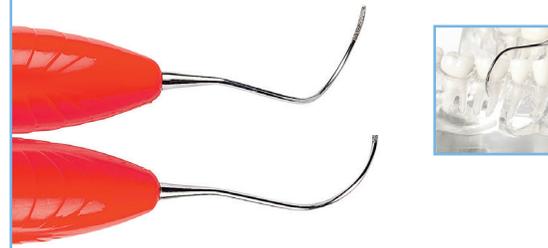

Furcation Diamond File

- Extremely delicate diamond coated instrument for definitive calculus removal in furcations and root depressions
- Easy access with a fine, flexible blade design

- LMErgoSense®** LM 261-262D ES
- LMErgoSense®** LM 261-262D EST
- LMErgoMax™** LM 261-262D XSI

Unique intelligence with **LM DTS™** RFID tag inside

To order add letter **T** to LMErgoSense® product code. E.g. LM 301-302 **ES** → **EST**



Handle options:



LMErgoSense®



LMErgoMax™



LMErgoNorm™

FURCATION INSTRUMENT

Well-suited for the cleaning of the furcations of multi rooted teeth.



Furcator KS

- Delicate, angled excavator
- Suitable also for the cleaning of narrow approximal concavities and grooves in incisors and premolars
- For cleaning of the furcation roof and inner surfaces, particularly concavations and grooves in them

- LMErgoSense®** LM 259-260 ES
- LMErgoSense®** LM 259-260 EST
- LMErgoMax™** LM 259-260 XSI

HOE SCALERS

For the removal of supra- and subgingival calculus.

Suitable for deep, narrow pockets and concave root surfaces. For both vertical and horizontal techniques.



Design:

- Fine hoe-shaped blade with rounded corners
- Elliptical cutting edge



Hoe Scaler, anterior

- For the anterior labial and lingual surfaces
- Works well also with lower incisors that are orally inclined



- **LMErgoSense®** LM 112-113 ES
- **LM 112-113 EST**
- **LMErgoMax™** LM 112-113 XSI
- **LMErgoNorm™** LM 112-113 SI




Hoe Scaler, lateral

- For the buccal and lingual surfaces of all teeth
- Can also be used for furcations



- **LMErgoSense®** LM 134-135 ES
- **LM 134-135 EST**
- **LMErgoMax™** LM 134-135 XSI
- **LMErgoNorm™** LM 134-135 SI

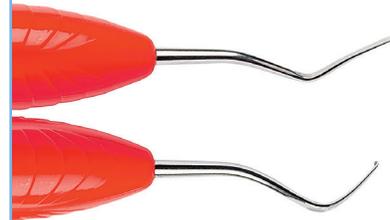

Hoe Scaler, posterior

- For the mesial and distal surfaces of molars
- Can also be used for furcations



- **LMErgoSense®** LM 156-157 ES
- **LM 156-157 EST**
- **LMErgoMax™** LM 156-157 XSI
- **LMErgoNorm™** LM 156-157 SI



SURGICAL HANDLE



LM-ErgoHold™ 3

- Excellent grip with elastic silicone coating
- Precise blade control
- Controlled rotation
- Blade position responsiveness maintained
- Controlled movement through tissue
- Fits all blades with fitment No. 3 (blades 10, 10A, 11, 11K, 11P, 12, 12D, 15, 15C and 15S)

- **LMErgoHold™** LM 9003
- **LM 9003T**

Unique intelligence with LM DTS™ RFID tag inside



To order add letter T to LM-ErgoSense® product code. E.g. LM 301-302 ES → EST



Handle options:



■ **LMErgoSense®**



■ **LMErgoMax™**



■ **LMErgoNorm™**